Now at your dealer or **wh.com**



Air scaler fast, gentle, powerful





Prophylaxis with style

Fast, gentle, powerful.

The air-driven scaler with the integrated spray and two large optic outlets is ideal for the removal of plaque and calculus as well as cleaning root surfaces. The result: your patients feel better in general and more self-confident thanks to their clean teeth surfaces.




Stepless regulation

The powerful, controlled deployment of the power can be adjusted to suit the application in question using the power regulator.



Large optic outlets

Two large optic outlets provide perfect illumination of the treatment site.

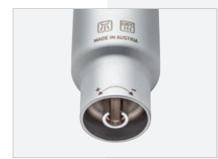

Tip changer

Simple changing of tips and protection against injury with the W&H tip changer.



Integrated spray

The integrated spray continuously cools the treatment site and simultaneously washes away the abraded microparticles and deposits.



First-class hygiene

Their smooth surface means that the instruments are easy to clean and can be both thermo washer disinfected and sterilized.



W&H adaptor set

The W&H adaptor set makes it possible for you to clean and disinfect the tips automatically using a validated procedure in the thermo washer disinfector – inside and outside.

Tips

for every application

Diverse and ergonomically perfect.

The diversity of the tip range and the optimal ergonomics of each individual tip make it easy for you to remove calculus gently and thoroughly even in hard-to-reach places.



Prophylaxis Tip 1AU

For removal of calculus, universally applicable, also for curved molar surfaces inclined towards lingual.



Prophylaxis Tip 2AU

For anterior teeth of the lower jaw, lingual surfaces and in cases of stubborn adhesion of thick, hard collars of dental calculus.



Prophylaxis Tip 3AU

For removing plaque from the neck of the tooth and removing nicotine stains.



Periodontology Tip 1AP

Removal of subgingival deposits, especially suitable for deep periodontal pockets.



Periodontology Tip 2APr. 2API

Tip for removal of subgingival deposits. Right-curved and leftcurved models available.



Periodontology Tip 3AP

With a diamond coating on the instrument tip for periodontal debridement of furcations and concavities.

Pictures are for illustrative purposes only. Additional equipment and accessories shown are not included as standard.

20403 AEN Rev. 003 / 31.01.2019 Subject to alterations

Technical data

Proxeo air scaler

| Туре: | ZA-55 L / ZA-55 LM / ZA-55 LS | ZA-55 / ZA-55 M | |
|---------------------------|--|---|--|
| Light: | yes | no | |
| Power regulation: | yes | no | |
| Oscillating frequency: | 6 kHz | | |
| Oscillation amplitude: | < 200 μm | | |
| Operating pressure: | $3~\pm0.3$ bar / $3~^{+0.6}_{-0.3}$ bar / $3~\pm0.3$ bar | $3~\pm0.3$ bar / $3~^{+0.6}_{-0.3}$ bar | |
| Coolant flow volume: | 30 – 60 ml/min | | |
| Coupling: | RQ-24, RQ-34 / RM-34 LED, Multiflex®* / Sirona®* | RQ-3, RQ-4, RQ-14 / Mulitflex®* | |
| Applicable W&H tips: | 1AU, 2AU, 3AU, 1AP, 2API, 2APr, 3AP Tips from other manufacturers must be approved for use with W&H air scalers. | | |
| In the scope of delivery: | Tip 1AU, 2AU, 3AU, 1AP, Tip changer | | |

 $^{^{\}star} \quad \text{Multiflex} \\ ^{\circledcirc} \text{ and Sirona} \\ ^{\circledcirc} \text{ are trademarks of third parties not affiliated with W\&H Dentalwerk B\"{u}rmoos GmbH.}$



| Type: | ZE-55 BC | ZE-55 RM |
|---------------------------|--|----------|
| Light: | no | |
| Power regulation: | no | |
| Oscillating frequency: | 6 kHz | |
| Oscillation amplitude: | <200 μm | |
| Operating pressure: | 2.2 – 2.5 bar | |
| Coolant flow volume: | min. 30 ml/min | |
| Coupling: | Fixed connection | |
| Applicable W&H tips: | 1AU, 2AU, 3AU, 1AP, 2API, 2APr, 3AP Tips from other manufacturers must be approved for use with W&H air scalers. | |
| In the scope of delivery: | 2 x Tip 1AU | |



W&H coupling RM-34 LED for Multiflex®* connection

- Use of LED without requiring changes to the dental unit
- Long LED service life
- Continuous regulation of the spray water supply
- Resuction stop
- 24-month warranty

Ignaz-Glaser-Straße 53, Postfach 1 5111 Bürmoos, **Austria t** +43 6274 6236-0 **f** +43 6274 6236-55 office@wh.com **wh.com**