

The new generation of rotary polishing







Tiny but mighty

The Proxeo Twist LatchShort Polishing System is a particularly smart solution for Prophy polishing: The small head with the optimally matched Prophy Cups significantly reduces the working height. You gain as much as four millimetres compared to standard latch systems.



Access all areas The molar region is more easily accessed from a distal to a mesial direction.



Comfortable and simple The easier the access the more comfortable the treatment.




Perfect

harmony






The smaller the head, the shorter the better the view of the treatment site.

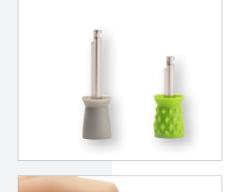


Measurably more space

The advantage of the new Proxeo Twist LatchShort system compared to the green contraangle handpiece: More space for polishing with the shorter shaft.

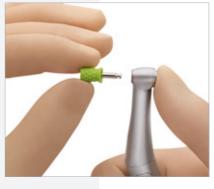

Change at the push of a button

With the practical push-button chuck system, you can change the Prophy Cups in seconds.

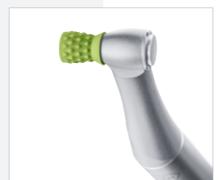

Perfect sealing

Saliva and polishing paste can no longer infiltrate the head so easily. This supports the long service life of the contra-angle handpiece.



the Cup and the thinner the neck,

Inspired by

Prophy professionals

In addition to the Proxeo Twist contraangle handpiece, W&H is now also offering Prophy Cups and Brushes. They were designed in consultation with Prophy users. Great emphasis was placed on features such as optimal adaptation, simple paste pickup and distribution as well as gentle cleaning up into the sulcus.

One-stop **shop**

The new Proxeo Twist LatchShort Prophy Cups in various hardness levels and Brushes are now available for optimal removal of biofilm and discolouration. All Proxeo Twist LatchShort Prophy Cups have a shorter shaft for an even lower working height.




Smooth surface at the end of the Cup for gentle polishing up into the sulcus.



special nub shape designed to prevent the paste from spraying.



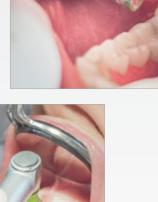

Nubs on the outside polish the interdental space at the same time.



Irresistible advantages for rotary polishing Slimmer. Finer. Innovative. The new Proxeo Twist LatchShort Polishing System with a working height that is up to 4 mm shorter thanks to the Proxeo Twist Prophy Cups and Brushes: for better access, more space and an optimal view of the treatment area.

Maximum protection and comfort



Slim and flexible

The contra-angle handpiece easily fits into the hand with its ergonomic design and slim neck. The combination with the W&H EM-12 L electric motor adds to its ergonomic properties.



Hygiene

The smooth exterior design makes the sterilizable and thermo washer disinfectable Proxeo Twist LatchShort contra-angle handpiece particularly easy to clean. The Prophy Cups and Brushes are disposable products and are not intended for reuse.

Technical data		
Туре		WP-66 W
Max. speed:		10,000 rpm
Transmission ratio:		4:1
Recommended speed for Prophy Cups and Prophy Brushes		1,250 – 2,500 rpm
Attachments that can be used:		W&H Proxeo Twist LatchShort Prophy Cups and Prophy Brushes
Prophy Cups and Brushes	(2)	Not for reuse

20480 AEN Rev. 002 / 28.01.2021 Subject to alterations



Melanie Ruckenthaler

W&H product manager

»The best of both worlds: the unique combination of a low working height and convenient push-button chuck system combines the advantages of screw and latch systems: easier access, better view and maximum comfort on the one hand and the ease of changing polishing Cups and Brushes on the other.«



Monika Riedl dental hygienist, Munich

»The Proxeo Twist LatchShort is the first system with a low working height and a push-button. With the working height as much as 4 mm lower than standard latch systems, there is more space and a perfect view of the treatment site.«



Implantology without Periodontology is like Yin without Yang.*



#NIWOP: The systematic and evidence-based workflow

Ignaz-Glaser-Straße 53, Postfach 1 5111 Bürmoos, **Austria** t +43 6274 6236-0 f +43 6274 6236-55 office@wh.com wh.com